

LIGO Laboratory / LIGO Scientific Collaboration

LIGO-T1000275-v1

LIGO

5/25/10

**Transmon SUS Telescope
Alignment Procedure**

Michael Smith

Distribution of this document:
LIGO Scientific Collaboration

This is an internal working note
of the LIGO Laboratory.

California Institute of Technology
LIGO Project – MS 18-34
1200 E. California Blvd.
Pasadena, CA 91125
Phone (626) 395-2129
Fax (626) 304-9834
E-mail: info@ligo.caltech.edu

Massachusetts Institute of Technology
LIGO Project – NW22-295
185 Albany St
Cambridge, MA 02139
Phone (617) 253-4824
Fax (617) 253-7014
E-mail: info@ligo.mit.edu

LIGO Hanford Observatory
P.O. Box 1970
Mail Stop S9-02
Richland WA 99352
Phone 509-372-8106
Fax 509-372-8137

LIGO Livingston Observatory
P.O. Box 940
Livingston, LA 70754
Phone 225-686-3100
Fax 225-686-7189

<http://www.ligo.caltech.edu/>

Table of Contents

1	Introduction.....	5
1.1	Off-axis Parabolic Design.....	6
1.2	Manufacturing Tolerance of the Off-axis Parabolic Mirrors.....	6
2	Alignment Procedure.....	6
2.1	Secondary Mirror Translation and Tilt Correction.....	6
2.2	Autocollimator and Shack-Hartmann Instruments.....	7
2.3	Initial Mechanical Alignment of Primary and Secondary Mirror Mounts.....	8
2.4	Coarse Alignment with Autocollimator.....	8
2.4.1	Initial Misalignment Conditions.....	8
2.4.2	Best Alignment with an Autocollimator.....	9
2.5	Alignment with the Shack-Hartmann Sensor.....	11
2.5.1	Commercially Available Shack-Hartmann Sensor.....	13

Figures

<i>Figure 1: ETM Transmon SUS Telescope.....</i>	<i>5</i>
<i>Figure 2: ZEMAX model of off-axis parabolic telescope.....</i>	<i>6</i>
<i>Figure 3: Double Pass Alignment Beam.....</i>	<i>7</i>
<i>Figure 4: Pick-off Beam Splitter for Returned Beam Sensing.....</i>	<i>7</i>
<i>Figure 5: Initial Misalignment Conditions.....</i>	<i>9</i>
<i>Figure 6: Wavefront after Alignment with the Autocollimator, 3-degrees of Freedom of Secondary Mirror.....</i>	<i>10</i>
<i>Figure 7: Spot Diagram after Alignment with Autocollimator.....</i>	<i>11</i>
<i>Figure 8: Alignment Using Five Degrees of Freedom of Secondary Mirror.....</i>	<i>12</i>
<i>Figure 9: Spot Diagram, after Alignment Using Five Degrees of Freedom of Secondary Mirror..</i>	<i>13</i>
<i>Figure 10: Wavefront Map from Shack-Hartmann Sensor.....</i>	<i>14</i>

Abstract

This technical note describes an alignment procedure for the Transmon SUS Telescope that uses a combination of an autocollimator for initial coarse alignment, and a Shack-Hartmann wavefront sensor for fine alignment. The telescope is part of the Transmon Suspension System for aLIGO.

1 Introduction

The purpose of this technical note is describe an alignment procedure for the ETM Telescope that uses a combination of an autocollimator for initial coarse alignment, and a Shack-Hartmann wavefront sensor for the final fine alignment. The telescope is part of the Transmon Suspension System for aLIGO.

A preliminary design layout of the Transmon Suspension with the folded ETM Telescope placed below the Transmon optical table is shown in Figure 1.

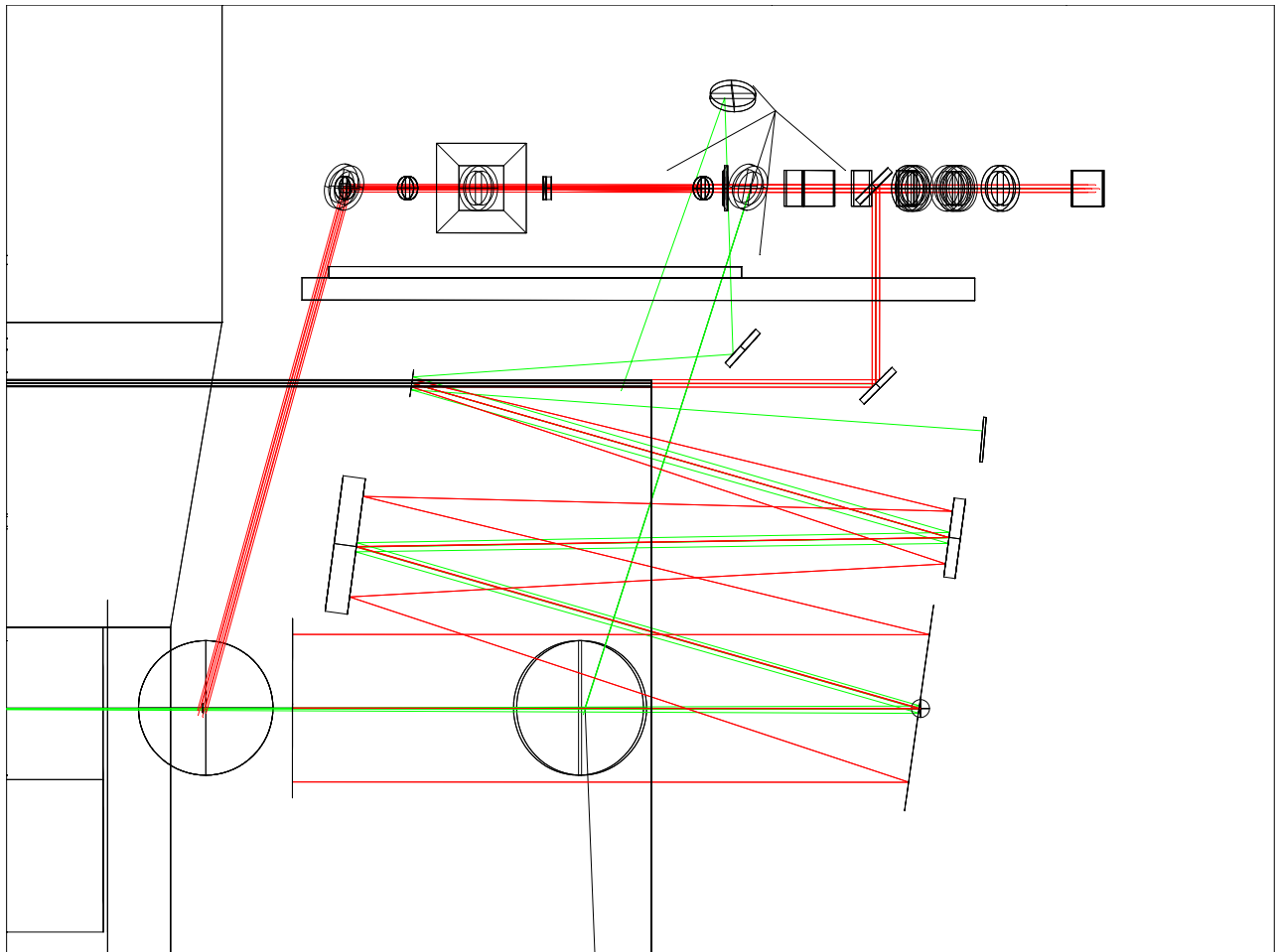


Figure 1: ETM Transmon SUS Telescope

1.1 Off-axis Parabolic Design

An unfolded ZEMAX layout was used to analyze the performance of the off-axis parabolic telescope, shown in Figure 2. The primary focal length is 2000 mm, and the secondary focal length is 100 mm. The off-axis distance for the primary is 562.2 mm, and for the secondary is 28.1 mm.

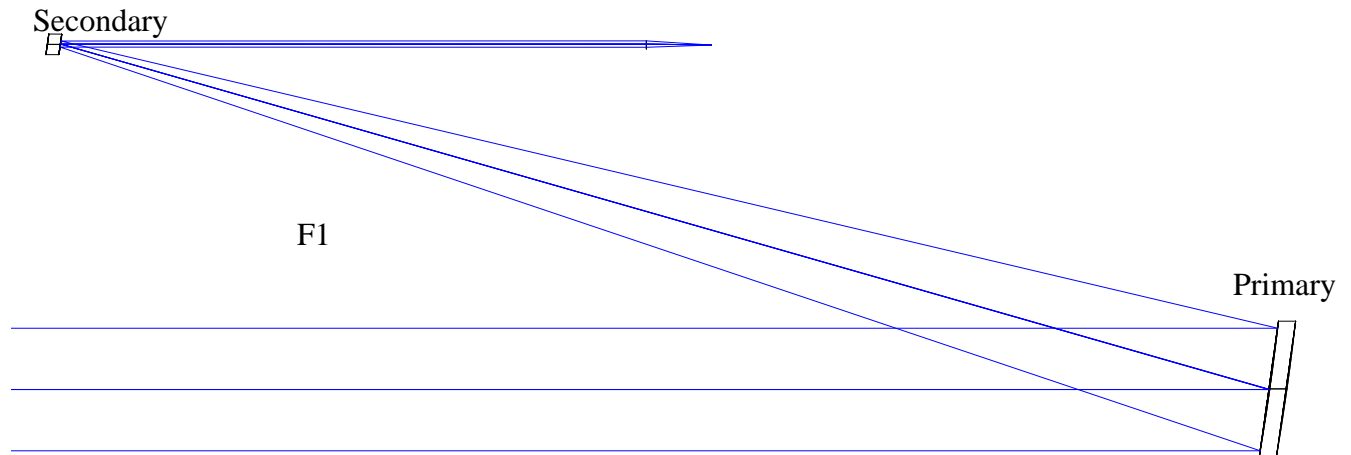


Figure 2: ZEMAX model of off-axis parabolic telescope

1.2 Manufacturing Tolerance of the Off-axis Parabolic Mirrors

The surface irregularity of the parabolas is specified to be $< 1/8$ wave @ 633 nm wavelength. This is equivalent to an irregularity of $< 1/13$ wave @ 1064 nm.

The tilt tolerance of the mechanical axis of the mirrors to the optical axis is ± 0.0083 deg.

The decentering tolerance of the primary mirror is ± 1 mm, and the decentering tolerance for the secondary mirror is ± 0.3 mm.

The focal length tolerance of the primary mirror is ± 20 mm.

2 Alignment Procedure

2.1 Secondary Mirror Translation and Tilt Correction

Satisfactory alignment of the telescope can be accomplished by moving only the secondary mirror in five degrees of freedom—tilting in pitch and yaw, decentering vertically and horizontally, and moving along the optical axis for focus.

The primary mirror will be placed in a fixed mirror mount whose mechanical axis is pre-aligned to within 0.1 mrad of the reference optical axis of the telescope. The rms combined error of the mirror mount and the tilt error of the primary is 0.01 deg.

The secondary mirror has the necessary degrees of tilt and translation freedom.

2.2 Autocollimator and Shack-Hartmann Instruments

Either an autocollimator beam, or a collimated light source, such as a laser or a super luminescent diode, will be injected into the secondary end of the telescope along the optical axis, as shown in Figure 3. The beam will reflect from a mirror when it emerges from the primary mirror and be reflected back through the telescope and emerge after a round trip at the secondary mirror. The beam splitter shown in the enlarged view will reflect a portion of the return beam to the cross-hairs of the autocollimator or to the Shack-Hartmann wavefront sensor.

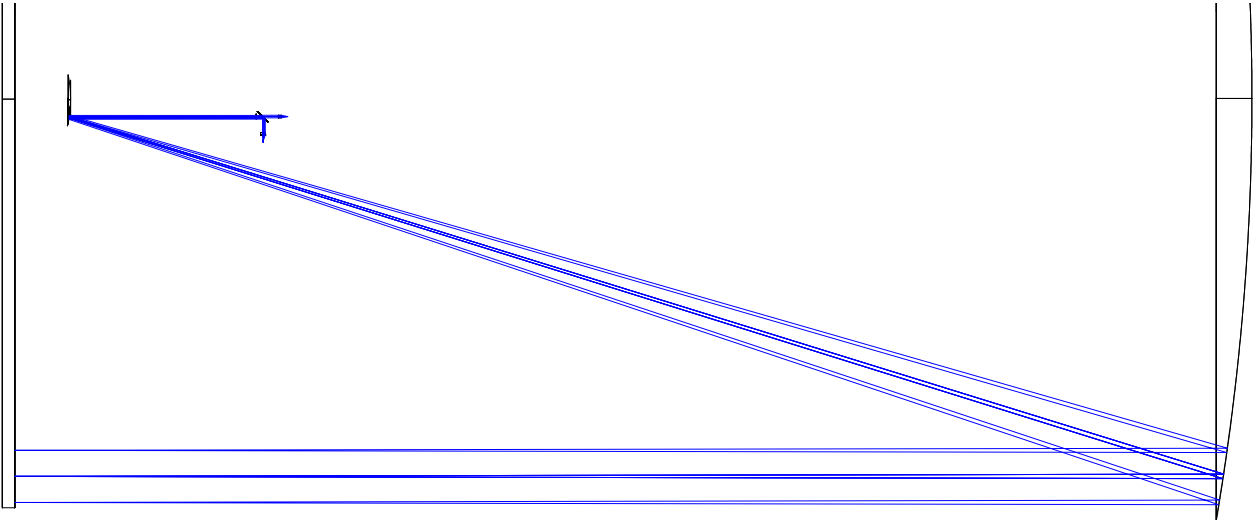


Figure 3: Double Pass Alignment Beam

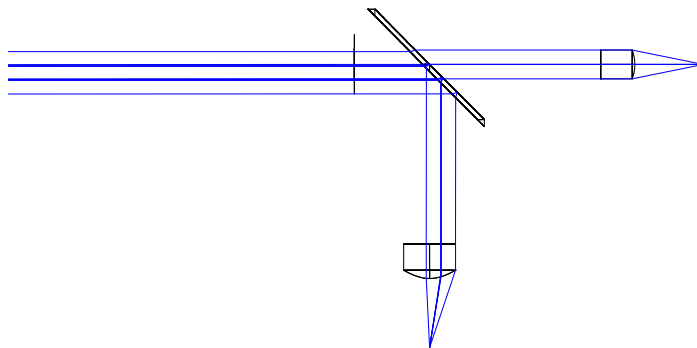


Figure 4: Pick-off Beam Splitter for Returned Beam Sensing

2.3 Initial Mechanical Alignment of Primary and Secondary Mirror Mounts

The ETM Telescope structure provides a reference surface whose normal defines the direction of the input and output optical axes of the telescope. The primary and secondary mirror mounts are pre-aligned parallel to the reference surface.

2.4 Coarse Alignment with Autocollimator

During this procedure, the secondary mirror will remain oriented in its initial mechanical alignment and will not be tilted.

1. Focus is achieved by moving the secondary mirror until the returned reticle pattern seen in the eyepiece of the autocollimator has sharp edges. After initial focus is achieved, it will be observed that the reticle pattern is not centered with the crosshairs of the autocollimator, which have been previously set coincident with the optical axis of the telescope.
2. Next, the returned autocollimator beam is aligned with the optical axis by translating the secondary mirror transverse to the beam until the reticle pattern is centered on the cross hairs.
3. The focus and beam pointing steps are iterated until the reticle pattern is visually in sharp focus and is centered on the cross hairs.

2.4.1 Initial Misalignment Conditions

The ZEMAX wavefront map in Figure 5 shows the initial aberrations of the wavefront of the telescope beam with the following worst case combination of tilt and de-center errors of the primary and secondary mirrors:

Primary

tilt	0.01 deg
decenter	1 mm

Secondary

tilt	0.01 deg
decenter	-0.3 mm

Defocus -20 mm

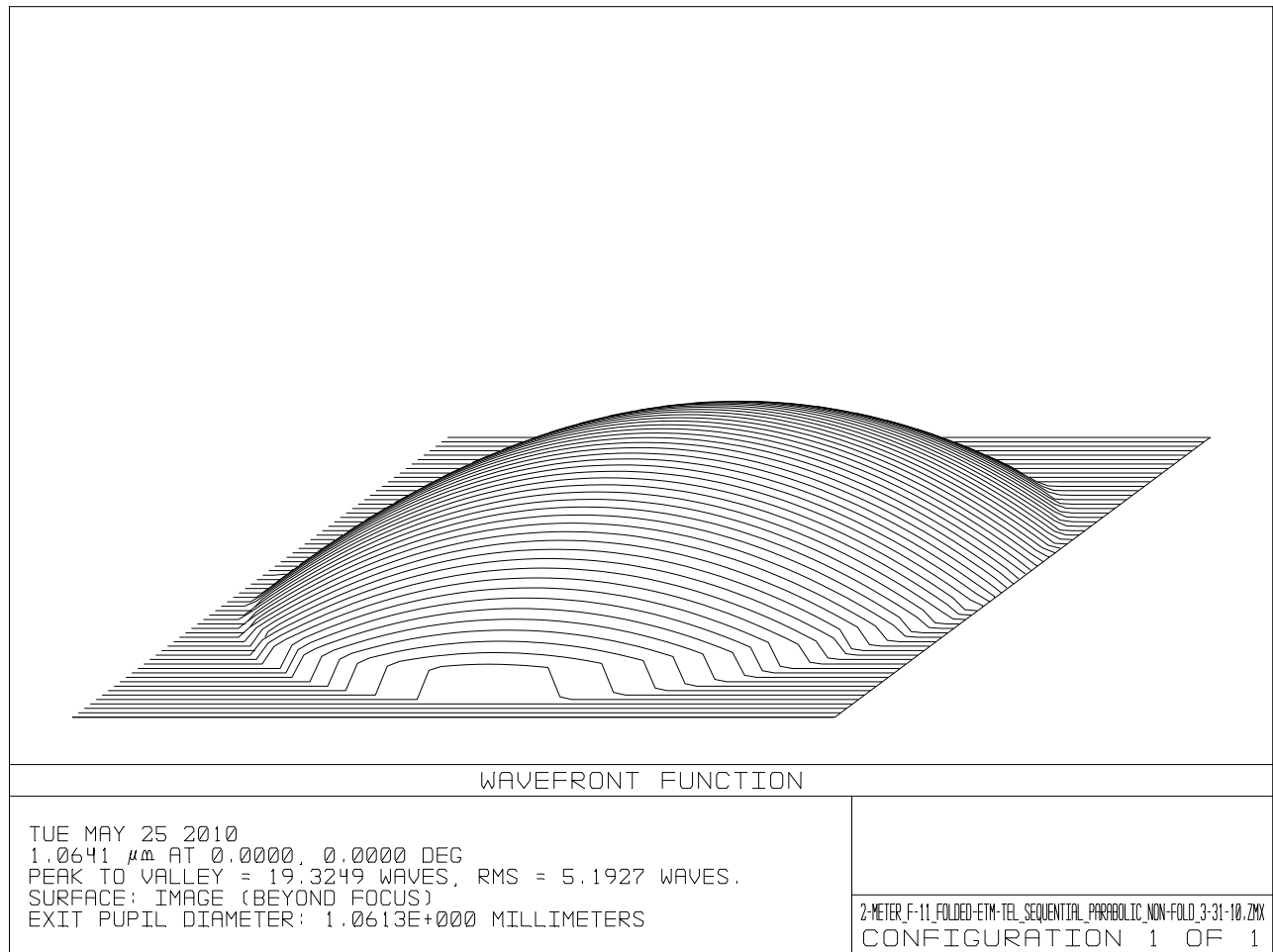


Figure 5: Initial Misalignment Conditions

2.4.2 Best Alignment with an Autocollimator

The resulting wavefront map after visually focusing and subsequently centering the reticle pattern on the cross hairs of the autocollimator is shown in Figure 6. However, it should be emphasized that this alignment was done by ZEMAX using a large number of trials seeking to minimize the wavefront distortion. The simple approach of focusing the reticle pattern by observing visually the sharpness of the edges will not achieve as good an alignment.

At best, this approach will result in a diffraction-limited the wavefront, but still exhibits a certain amount of astigmatism as seen by the saddle shape of the wavefront map in Figure 6. The astigmatic beam shape is exposed in the spot diagram of Figure 7. The black circles represent the Airy's disc diameter, and are approximately located 1 Rayleigh range before and after the focused spot. The spot diagram implies that the astigmatic foci are within the Rayleigh range of the telescope secondary beam, and may be problematic for the operation of the Transmon SUS system.

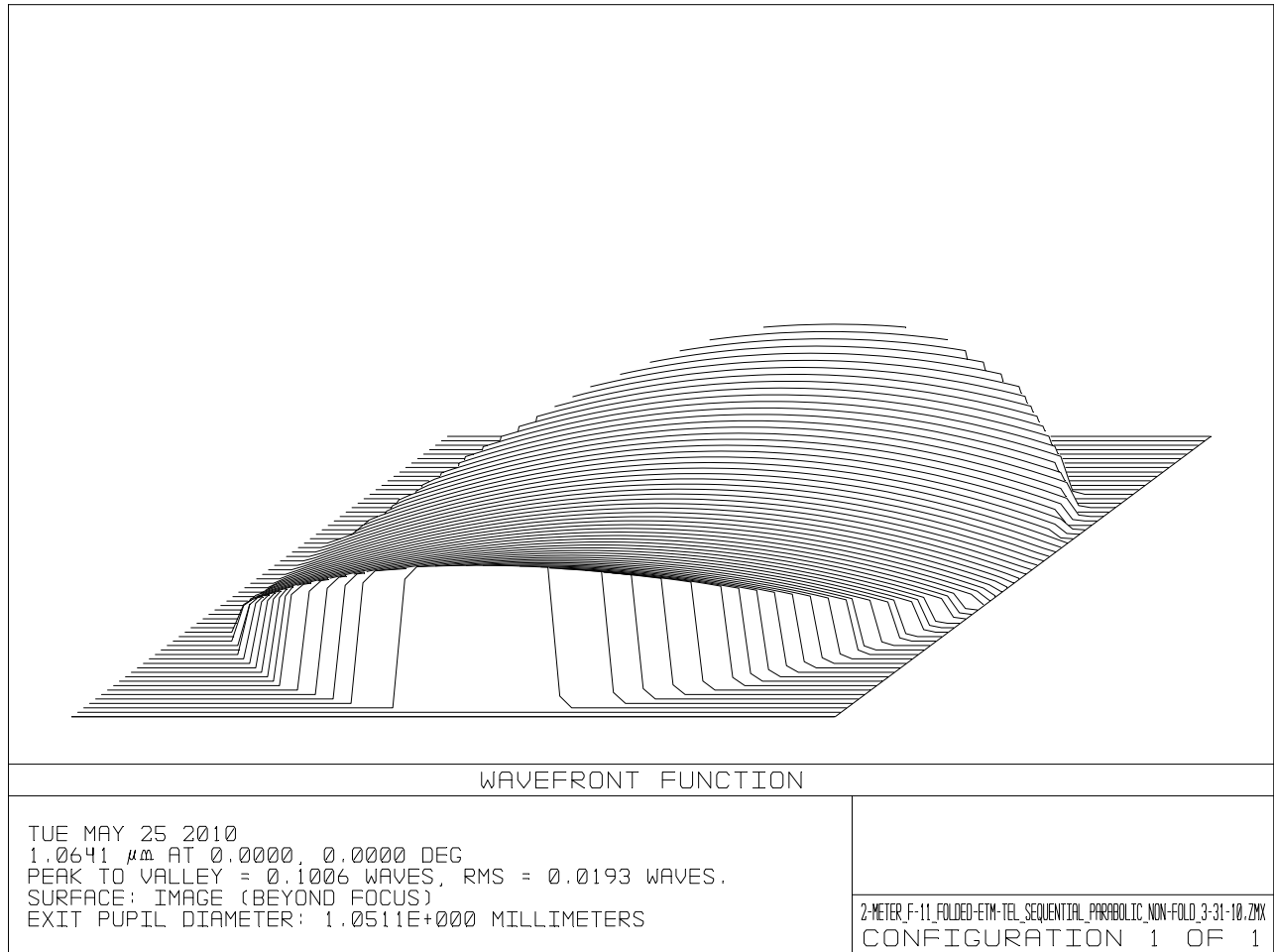


Figure 6: Wavefront after Alignment with the Autocollimator, 3-degrees of Freedom of Secondary Mirror

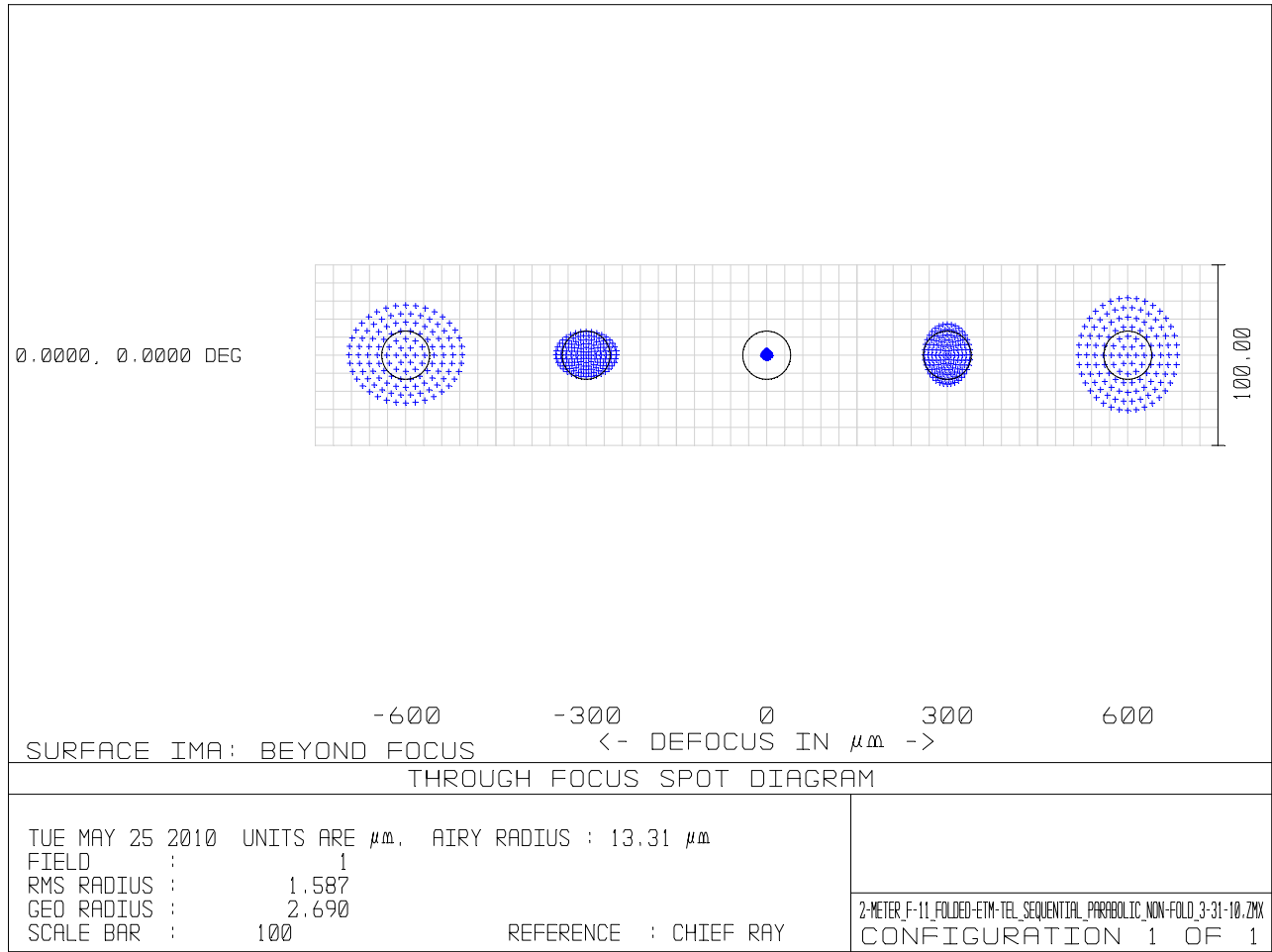


Figure 7: Spot Diagram after Alignment with Autocollimator

2.5 Alignment with the Shack-Hartmann Sensor

The wavefront distortion and the resulting astigmatism can be minimized by allowing the additional degrees of freedom of tilting the secondary mirror. The Shack-Hartmann sensor is needed to provide a direct measurement of the wavefront map, and to provide a quantitative tool for aligning the five degrees of freedom of the secondary mirror position, which would be impossible to do with the autocollimator.

The resulting ZEMAX alignment using five degrees of freedom is shown in Figure 8. The wavefront aberration is < 0.004 waves, which is less than the specified $1/13$ waves wavefront irregularity for the individual telescope mirrors.

The spot diagram in Figure 9 also shows negligible astigmatism.

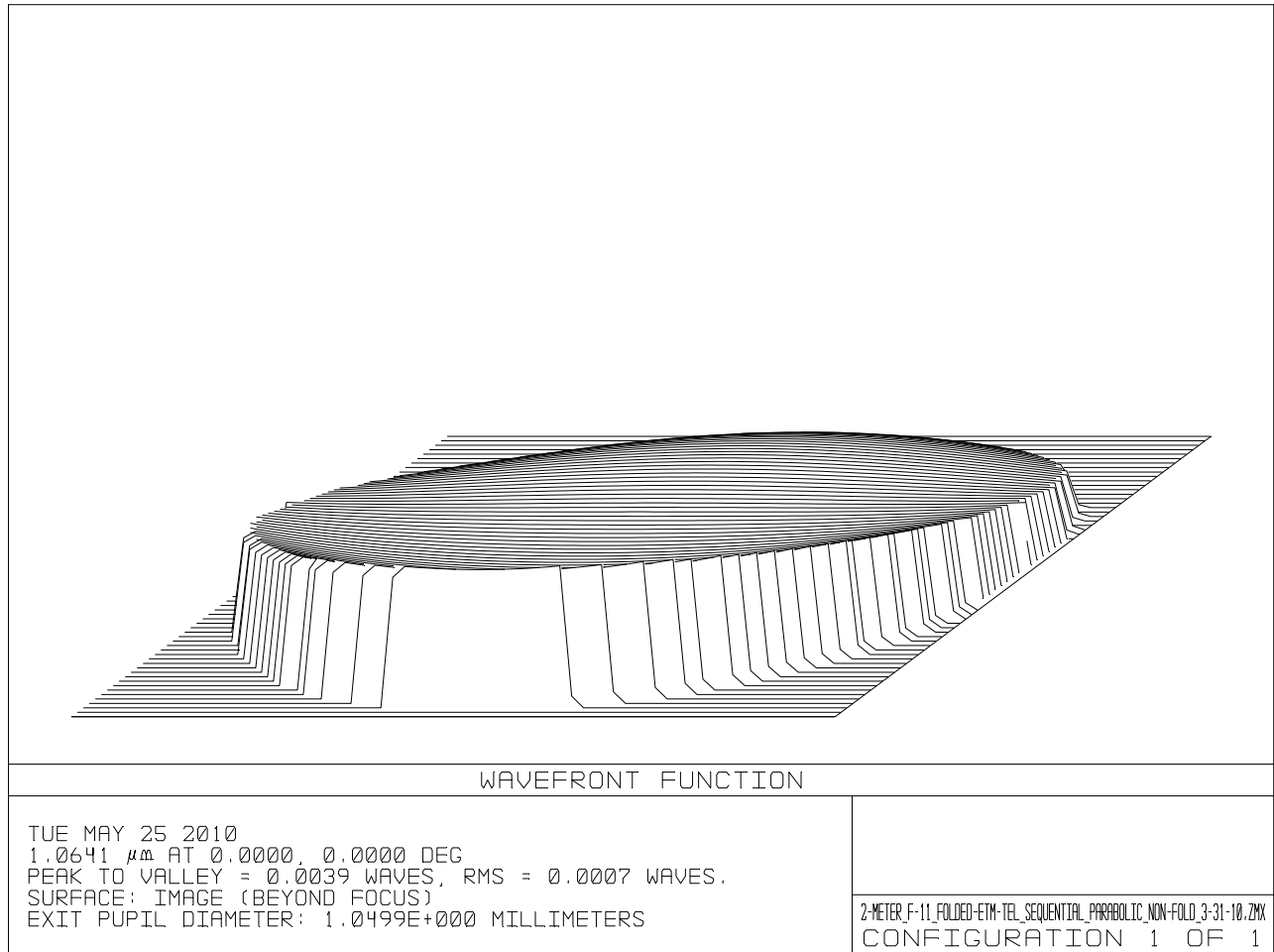


Figure 8: Alignment Using Five Degrees of Freedom of Secondary Mirror

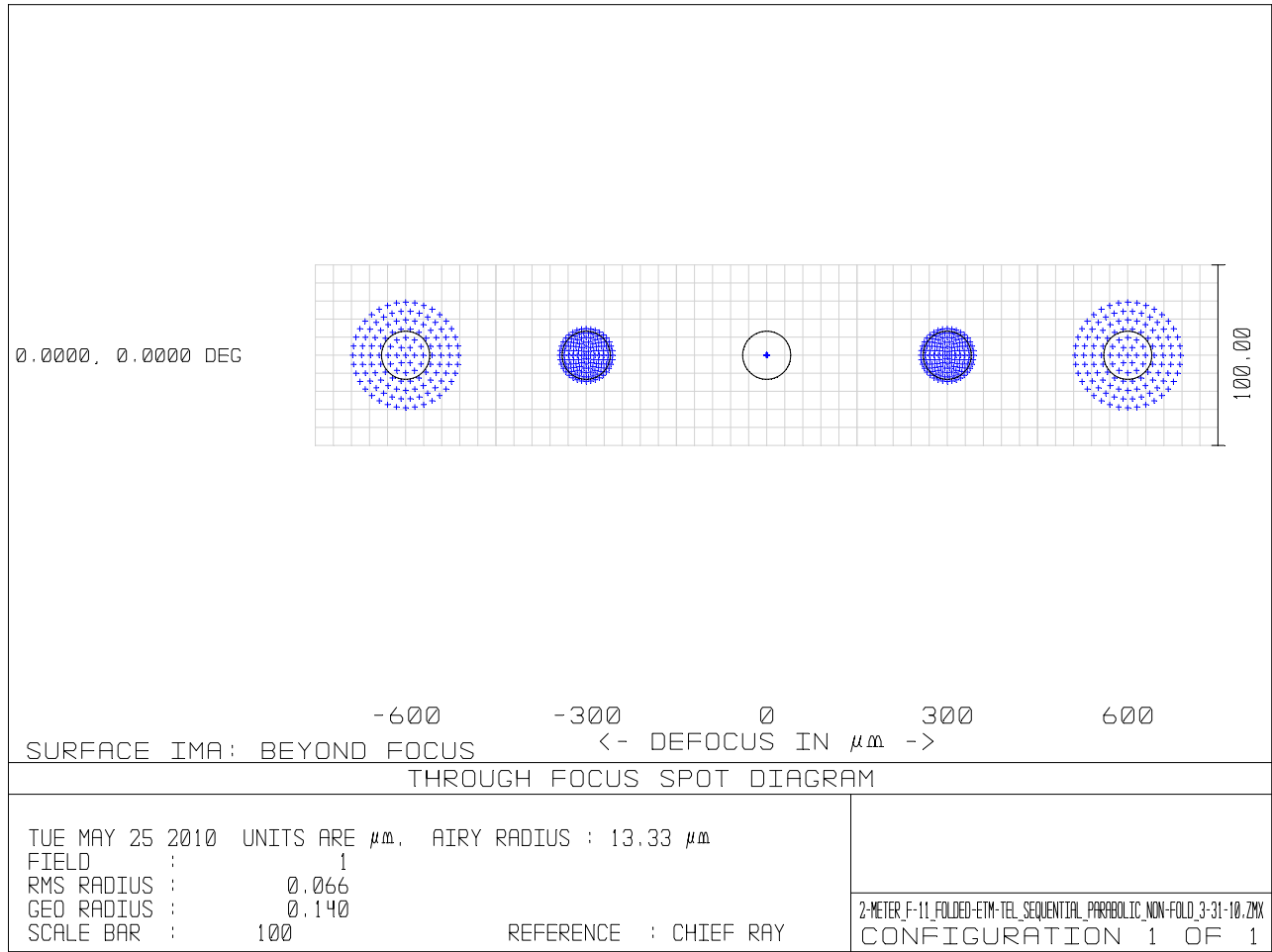


Figure 9: Spot Diagram, after Alignment Using Five Degrees of Freedom of Secondary Mirror

2.5.1 Commercially Available Shack-Hartmann Sensor

A screen capture of the output of a commercial Shack-Hartmann sensor is shown in Figure 10. It provides quantitative data of the wavefront distortion in terms of Zernike coefficients.

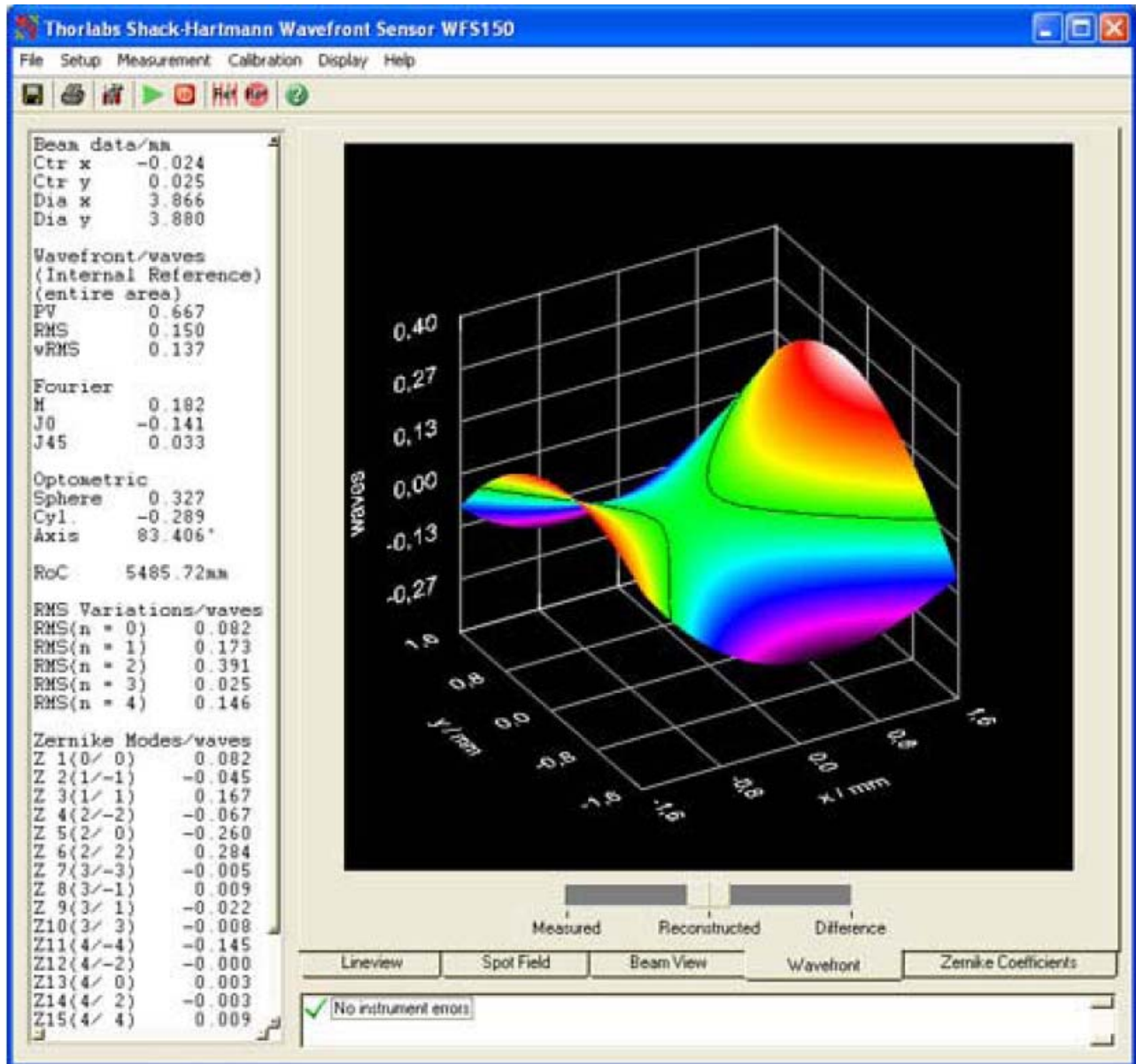


Figure 10: Wavefront Map from Shack-Hartmann Sensor